

VACUUM INTEGRATED PUMP

VIP®



DESCRIPTION

The ALPS VIP® is an air expulsion integrated pump for elevated vacuum for prosthetic socket. The use of VIP provides patients continuous total contact when donning the prosthesis. This is essential for an even weight distribution and to distribute the suspension of the prosthesis over the whole surface of the residual limb. Elevated vacuum created by the VIP® stabilizes the residual limb, overcoming all issues connected with the residual limb changing volume during the day.

PART

VIP-DH

VIP-DL

TOP FEATURES & BENEFITS

-  Efficiently deals with stump volume fluctuation which in turn leads to improved prosthesis control and comfort
-  Reduces the intensity of peaks in pressure on the patient's limbs, improving comfort and providing greater freedom of movement
-  Helps preserve residual limb skin health by improving circulation and consequently increasing limb oxygenation
-  Does not affect the cosmetic appearance of the prosthesis

Easy to
install

Improves
blood flow

Improves
comfort

Improves
walking

Quiet

Light

AVAILABLE IN 2 DIFFERENT VACUUM LEVELS

CHOOSE LEVEL BASED ON:




ACTIVITY LEVEL




SKIN CONDITION

HIGH VACUUM (Red)
10-13 in of mercury

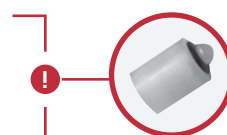


We highly recommend
using the **mini valve** to
prevent distal distraction


LOW VACUUM (Green)
5-8 in of mercury

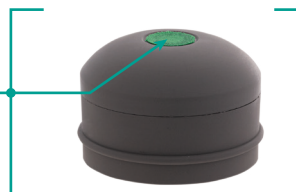

WITHOUT MOUNTING PLATE

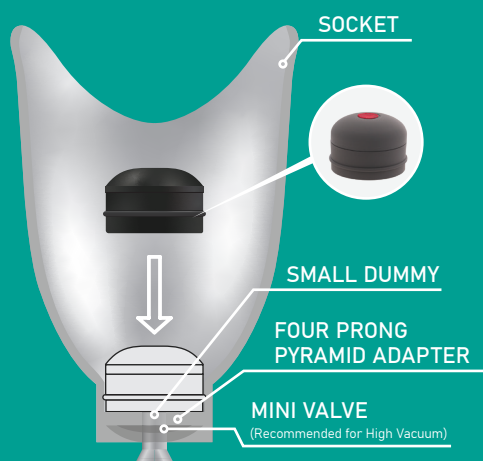

WITH MOUNTING PLATE

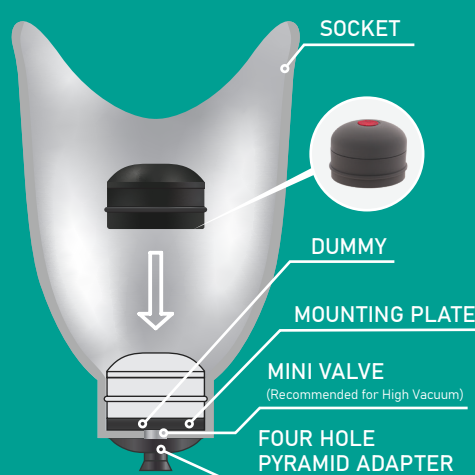

PRODUCT SPECS

	HIGH VACUUM (RED)	LOW VACUUM (GREEN)
PART #	VIP - DH	VIP - DL
DIAMETER	66MM	66MM
HEIGHT	48MM	48MM
WEIGHT	75G	64G
MATERIAL	SILICONE	SILICONE

USE CODES BELOW WHEN ORDERING

PART NAME	PART #
DUMMY	VE319 - D
SMALL DUMMY	VE319 - SD
MOUNTING PLATE	VE819 - P
LOCK (SMALL)	VIP-L1120-S
MINI VALVE (Highly Recommended)	MV-009
MINI VALVE KIT	LMV-009
PIN (4.8cm)	VIP-92L
PIN (3.0cm)	VIP-92S



Tel +420 377 223 127 | Fax +420 377 223 010

www.easyliner.eu | info@easyliner.eu

Božkovské náměstí 17/21, 32600 Plzeň – Czech Republic

