

PROXIMAL SEAL INTERFACE





PSI

Gel Liner, Distal Fabric – Without
Distal Attachment

DESCRIPTION

ALPS™ PSI Liner is formulated with EasyGel and provides comfort, strength, and elasticity for a variety of amputees. The PSI can be used in conjunction with most valve suction or vacuum systems used in the market today. The fabric cover, located distally, allows the air to wick out of the prosthesis. This allows a sleeve to seal on the liner just proximal to the trim lines of a transtibial socket.

TOP FEATURES

-  Exposed gel portion allows for means of sealing
-  Fabric located distally to act as air wick
-  Unlimited restriction of stretch of the proximal portion of the liner
-  Available in 6mm uniform and tapered from 6mm at distal end to 3mm



EASYGEL



CUSHION

MAKING LIVES BETTER




#makinglivesbetter




CROSS SECTION DIAGRAM

6.0MM UNIFORM OR 6/3MM TAPERED

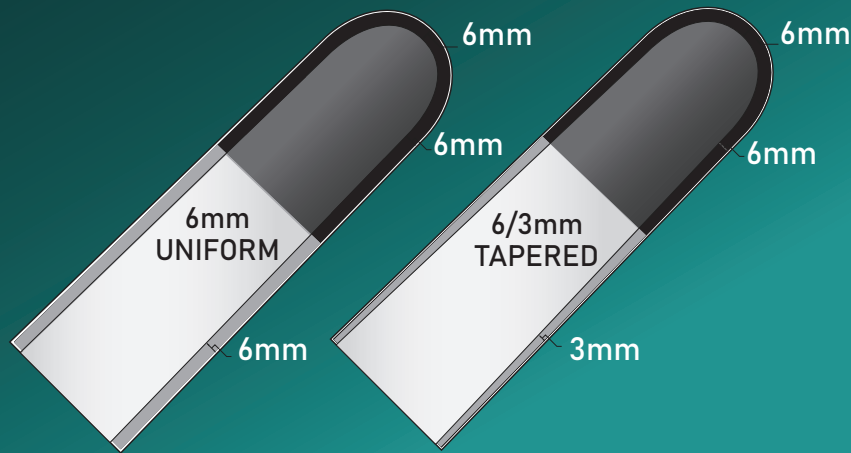

HOW TO SIZE

Measurements taken
6cm from the distal end

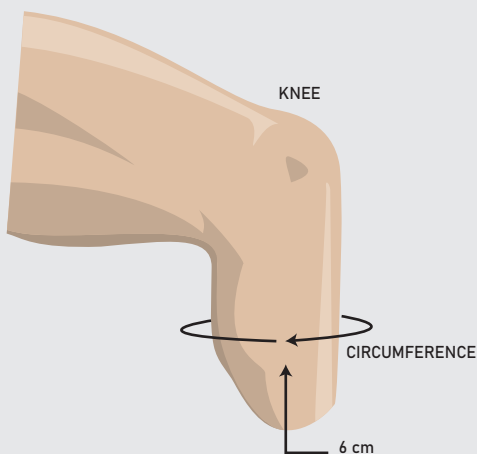

SIZING CHART

MEASURED SIZE (CM)	CUSHION 6MM UNIFORM OR 6/3MM TAPERED
16-19 cm	PSI16 - (6 or 6/3)
20-23 cm	PSI20 - (6 or 6/3)
24-25 cm	PSI24 - (6 or 6/3)
26-27 cm	PSI26 - (6 or 6/3)
28-31 cm	PSI28 - (6 or 6/3)
32-37 cm	PSI32 - (6 or 6/3)
38-43 cm	PSI38 - (6 or 6/3)
44-53 cm	PSI44 - (6 or 6/3)



Tel +420 377 223 127 | Fax +420 377 223 010

www.easyliner.eu | info@easyliner.eu

Božkovské náměstí 17/21, 32600 Plzeň - Czech Republic

MAKING
LIVES
BETTER